



NUA SKIN TAG REMOVAL SYSTEM SKIN TAG/MOLE REMOVAL PEN INSTRUCTIONS

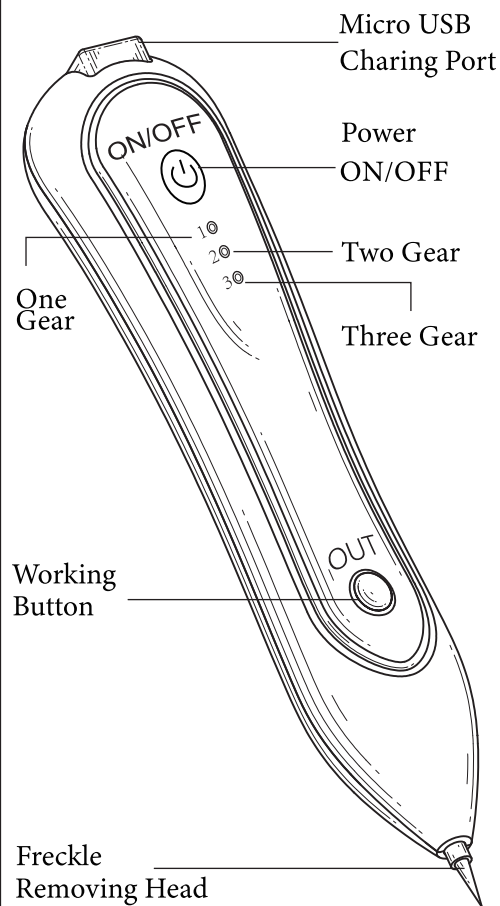



Mole Removal pen

Featuring carbon-ionization technology, this Skin Tag Removal Pen can safely, effectively, and bloodlessly remove moles and freckles without receiving any uncomfortable electric shock.

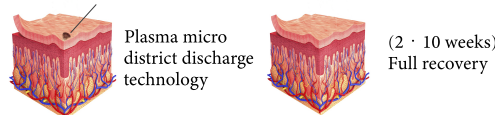

Needle Comparison:

1. Coarse needle (for repeated use): Use on larger moles, freckles, warts, etc. This more robust and powerful needle can remove blemishes quickly.
2. Fine needle (one-time use): For use on relatively small moles, freckles, and other blemishes. The fine needle is less painful and is better for detail work.

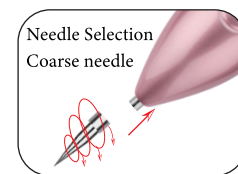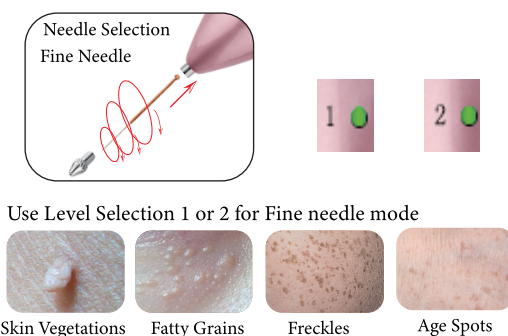

Use Level Selection 3 for coarse needle mode



Precautions: Before Use

1. Before use, disinfect the needle and the target area of skin with an alcohol pad.
2. For sensitive individuals or when using higher powers, numbing cream (not included) can be applied to the target area before removing blemishes.
3. Before removing your first mole or freckle, it is advised to practice for a short time on a less sensitive area in order to familiarize yourself with the pens' functions, as well as with the sensation during use.
4. Determine Strength Level:
 - A. For smaller freckles and moles, use setting level 1 and/or 2.
 - B. For larger blemishes, use setting level 3.
5. For proper use, only touch this pen lightly to each spot, targeting only the top layer of the skin.
6. It is recommended to gently sweep back and forth across the blemish.
7. For deeper moles, tilt the pen and use a lower strength with multiple applications.
8. Do not pierce the skin with the needle.
9. Do not use on one area/blemish for too long.

Precautions: After use

1. Patience is crucial.
2. One treatment /use at a time. Wait until treated area has healed from first treatment before applying a second treatment to the same area.
3. During the healing period, keep exposure to the sun at a minimum.
4. For small to medium sized blemishes after one week, a scab will have formed over the treated area. After about 2 weeks, the scab should come off naturally, and all trace of the blemish will be gone.
5. For larger moles and freckles, healing may take longer. Any remaining discoloration means that the pigment is deeper in the skin. In this case, the area can be treated again, but make sure to wait until the first treatment has healed completely.
6. In addition, taking care of your general health by staying hydrated and consuming foods high in healthy oils will contribute to a speedy recovery. I.E (drinking plenty of water, eating more vegetables, citrus fruits and nuts)
7. On the other hand consuming foods such as ginger, beef, and foods high in salt such as soy sauce is not recommended during the recovery period.
8. Don't drink wine or take aspirin (or similar drugs) within 1 week of a blemish treatment. As these contain chemicals that thin your blood.
9. Strenuous exercise that causes excessive sweating should be avoided within 1 week of mole/freckle removal to avoid causing skin to scar.
10. Application of some sort of healing cream such as cortisone or aloe Vera can further aid recovery.

Precautions

1. For first-time users, it is recommended to test the device on the back of your hand before using on a sensitive, delicate area such as your face.
2. Keep out of the reach of children at all times.
3. Avoid water on the wound for 1-2 weeks. If wound remains too humid it may become inflamed.
4. Don't pick off the scabs left by the blemish treatment.
5. Be patient and let nature take its course.
6. Depending on the size of the blemish, recovery periods can range from 2-10 weeks.

Usage Instructions

Remove Spot:

For removing spots, please place the mole removal touching the spot slightly, gently sweep back and forth, it can be clearly found that melanin will burn and fall off.

For small spots, it can be removed only one time.

For the darker spot, it's suggested that carbonizing be done 2-3 times. Remain as still as possible, and keep hand holding the device balanced during the process, otherwise it will cause burning, it also affects the wound recovery, resulting in scar.

It is recommended that the person performing the removal BE OTHER THAN the person having the blemish removed.

Remove Dot mole:

When removing dot mole, please pay attention to distinguish the extent of moles. Some moles grow on skin epidermis, while others are shallow. Do not pierce the surface of the mole at once, it's better to burn in small batches.

For the bad mole, please avoid inflammation of the wound after treatment, it's preferred to put some numb cream for the first 2-3 days of the wound recovery.

Note: Stop using numb cream if the wound becomes inflamed.

Remove Granulation:

Please lightly touch the granulation with pen tip to burn the derived granulation cells. They will automatically fall off over a period of one week.

For the small granulation, it can be removed only one time, while multiple operations to remove the larger or tougher granulation in larger areas should be used to avoid longer recovery.


Remove Tattoos:

Note: This device is only recommended to be used on small tattoos, for tattoos larger than 5cm it is recommended to seek professional medical or cosmetologist intervention.


For small tattoos, multiple operations will be needed and this process should be done over a longer period of time, allowing for recovery periods of the area between uses. To avoid any scars and to avoid prolonging the recovery period between uses make sure to follow all precautions listed on this manual. For areas larger than 3cm, use the coarse needle.

For severe inflammation, please immediately treat in hospital.


Operation Steps:

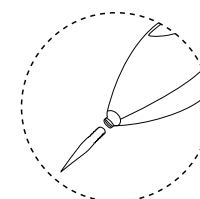


1. Choose the appropriate needle and install it on the pen point.

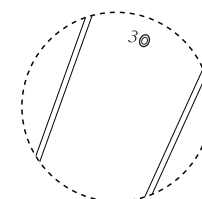

2. Long press ON/OFF button for 3 seconds to turn on the pen.

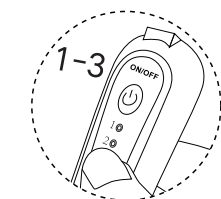

3. Click the ON/OFF button to select speed level.

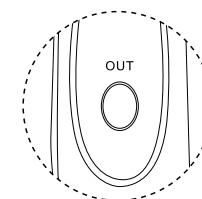

4. Press and hold OUT key to start working.

Manufactured For:
Home Aide Diagnostics, Inc.
1072 S. Powerline Rd.
Deerfield Beach, FL 33442

Distributed by:
Nu Skin Technology

For any questions please call: 1-800-915-0116

Made in China