

nuoskin FACIAL SCRUBBER

ULTRASONIC SKIN SCRAPER USER MANUAL



Thank you for purchasing this ultrasonic skin scrubber. To use this product correctly, please read this user manual before using it.




Product Structures & Parameters



Name	NUASKIN FACIAL SCRUBBER ULTRASONIC SKIN SCRAPER
Technology	Ultrasonic Ion Vibration Technology
Working Voltage	3.7V
Charge Current	5V/1A
Rated Power	3W
Mode 1	Skin Cleansing Mode
Mode 2	Nutrition Import Mode


Skin Cleansing Mode:

1. Apply a cleansing lotion or clean water evenly to the face for fully moisturize the face;
2. Press the switch to open the scrubber, and slowly move the scrubber along the facial skin. The dirt and grease on the skin are removed by the high-frequency ultrasonic wave.
3. After use, turn off the machine and clean the scrubber with a clean cotton cloth.

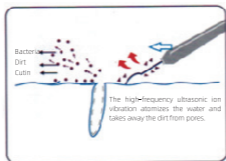

Import Nutrition Mode:

1. Apply an appropriate amount of water emulsion or serum to the face and neck;
2. Reverse the skin scrubber and open the switch. The reversed probe contacts the skin at an oblique Angle of 45 degrees, slowly move the probe to import nutrition into skin;
3. After skin care, turn off the skin scrubber and clean the scrubber machine with a clean cotton cloth.

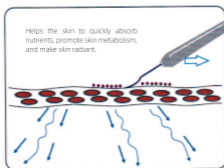

Product Usage Notes:

1. When using, the probe should be used at an angle of about 45° to the skin surface;
2. Please do not use this product in the throat, eyeball, mucous membranes, etc.;
3. Do not stop at some point in the face all the time when using, you need to use the probe to slowly move on the face.
4. Please clean the probe with a moist cotton pad before use.
5. Before using for the first time, please fully charge the machine. When Charging, the charging indicator will light red; when fully charged, the charging indicator will light green; it will take about 3 hours to fully charge.
6. The instrument has a set usage time. It will automatically turn off after about 10 minutes of use. If you need to continue using it, please restart the instrument.

Safety Precautions:

1. Please turn off the machine before cleaning the product;
2. This product is not suitable for use in humid environments. The probe can be washed with water, but the body part of the product is not waterproof. Do not immerse the whole machine in water;
3. Since the product has a built-in rechargeable battery, do not put the product in a fire or use it in a high-temperature environment;
4. If you are allergic to metal, please do not use this product;
5. When using, if the product is abnormally hot, please stop using it immediately and shut down to avoid burns;

Do not use this product in the following cases:

- A. During pregnancy, postpartum lactation futures women;
- B. People with allergies, skin diseases, skin inflammation, skin allergies, other skin diseases, and severe sunburn;
- C. People with skin damage, eczema, swelling, etc.;
- D. People who have just finished surgery and are recovering;
- E. When the body is unwell or has a fever;
- F. When drinking alcohol or taking drugs.