

# USER MANUAL

## NUASKIN VACUUM PRO BLACKHEAD REMOVAL DEVICE

---



Made in China


# Product Information

This multi-functional cleansing blackhead instrument uses a non-invasive vacuum negative pressure technology. 3 different suction modes are designed for 3 different skin types, which can deeply absorb blackheads, excess oil, acne, and dirt on the skin. To meet various usage needs, it's equipped with 6 different tips. It cleans the skin in all directions and gently lifts the skin, making the pores smaller, improving the secretion of skin oil, increasing the elasticity of the skin, and making the skin cleaner and smoother.

## Features:

- 3 Suction Modes:** According to different skin types, suitably designed suction modes for dry skin, oily skin, and mixed skin.
  - 6 Functional Suckers:** In order to meet the different needs of sucking blackheads, cleaning pore and microdermabrasion, It is specially equipped with 6 suckers with different functions.
  - High-definition LED Display Screen:** High-definition LED screen shows the working status of the product instantly, helping you to use the product easily.
  - Rechargeable Battery:** Equipped with high -capacity 380mA polymer batteries and USB charging cable. After fully charged, it can last for 90 minutes.
  - Low Noise:** The maximum suction level at work produces no more than 30 decibels of sound standard recommended by the International Organization for Standardization (ISO).
-

# Suitable For All Skin Types

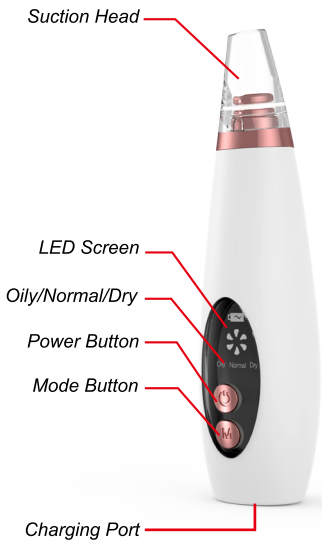

## 6 Replacement Heads



### Oval Head

Used to reduce wrinkles and increase skin elasticity



### Mini Round Head

Suitable for sensitive skin and less blackhead skin



### Round Head

Suitable for all types of skin



### Small Round Head

Suitable for sensitive skin and weak parts of the skin



### Middle Round Head

Used for oily skin, suck out blackheads and acne



### Big Round Head

Used to suck out blackheads and lift the skin

---

# Using Steps

1. **Charge:** Make sure to charge the device for 6 hours before use. (Red light indicates charging, green light indicates fully charged)
2. **Clean:** Make sure the skin is clean and dry before using this blackhead remover
3. **Before Use:** Use a hot towel and place it on your face for 3-5 minutes to open up the skin pores. You can also use a facial steamer or use the blackhead remover immediately after bathing.
4. **Using:** According to your skin, choose the appropriate probe and install it. Press the power button and choose the power lever (Oily/Normal/Dry) of suction force. Use in sweeping motion 1-2 seconds long.
5. **After Use:** Wash your face with cold water immediately after use. After washing face, you can put on a cold pack or rub an ice cube on the face to tighten pores.
6. **Clean Probe:** Properly clean the probe every time after use to sanitize and to keep from blocking the filter. You can clean the probe with 70% or less alcohol or wash with soap and a wet washcloth.

## Please Note:

\*This product may not be suitable for thin, sensitive or spider veins skin, please test on forearm before first use. If allergic or sensitive to the vacuum technology, please don't use it.

\*When using the product, please move its suction head in circular or linear direction slowly on your skin, and make sure it stays no more than 2 seconds for each skin area. Otherwise, it may hurt your skin and make it purple (Especially for sensitive skin) because the suction power is strong.

---

## FAQ

**Question 1:** Can I use the device when charging?

**Answer:** No, please do not use while charging. It's a safety precaution.

**Question 2:** Why does my face feel sunburned or windburned?

**Answer:** This device may cause redness and a feeling like a windburn, but it will go away after a few minutes. Do not worry, this is a normal situation. Wash your face with cold water and apply essence and moisturizer. You can also follow up with rubbing an ice cube all over the face to close pores and calm skin.

If you experience bruising or redness lasting hours then it was because you are keeping the vacuum probe too long in one place please refer to the right using steps and make sure it stays no more than 2 seconds for each skin area.

**Question 3:** Can I use the device when wet skin?

**Answer:** Please do not use on wet skin. The wet skin does not allow for a proper seal to be formed by the vacuum and will not provide optimal results. All electrical equipment must be used in a dry environment to avoid safety hazards.

**Question 4:** Can I use the device and blackhead remover gel and kit tool together?

**Answer:** Yes, before using the product, use a hot steamer or hot towel or blackhead remover gel/cream on your skin for 5 minutes to completely open your pores first. For some hard-to-remove blackhead, we suggest you use blackhead remover gel and kit tool together for best performance.

---

# Product Parameter

Product Name	Blackhead Remover Device	Working Voltage	5V
Suction Strength	$\leq -62\text{KPa}$	Charging Voltage	5V
Using Time(Fully Charged)	90 minutes	Input Power	2.5W
Certification	CE & ROHS	Rated Voltage	3.7V
Battery Type	Polymer Battery	Battery Capacity	380 mA

# Packaging Accessories

1



Remove Blackhead Device

2



Suction Heads\*6

3



USB Cable

4



User Manual

